
User manual

Gazelle PS

R82[®]
by Etac

M1110 - 2024-06-07 - Rev. 14.0

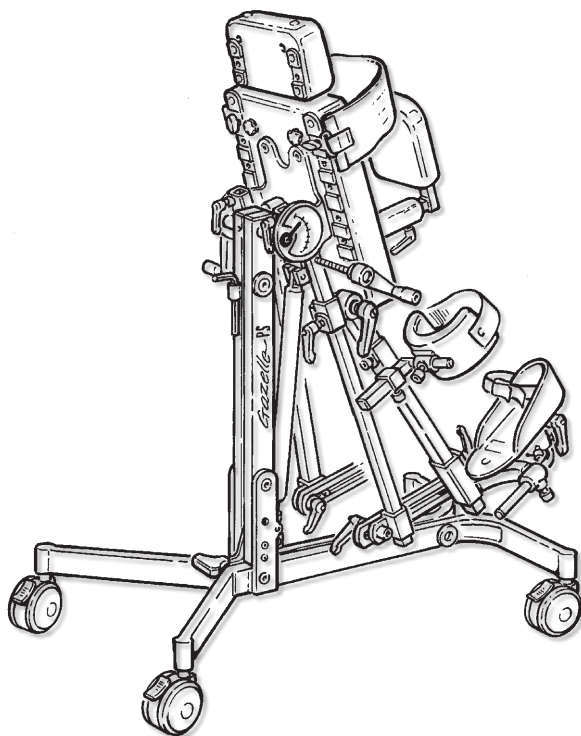www.etac.com


	3	ENGLISH DEUTSCH	12-21 22-33
	4	NEDERLANDS DANSK	34-43 44-53
	5	NORSK SVENSKA	54-63 64-73
mm (inch) kg (lb)	6	SUOMI ESPAÑOL	74-83 84-93
	7	FRANÇAIS ITALIANO	94-103 104-113
	8 - 11	PORTUGUÊS ČESKY	114-123 124-133


1.


Gazelle PS 1

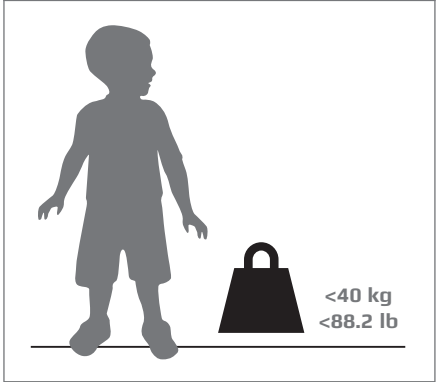
Gazelle PS 2

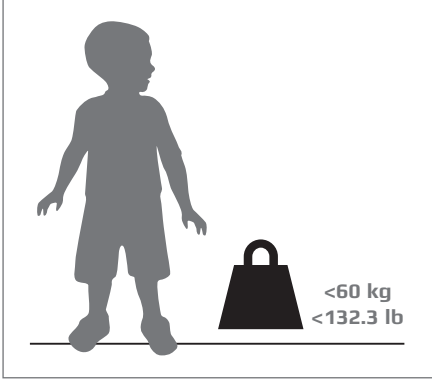Gazelle PS 3

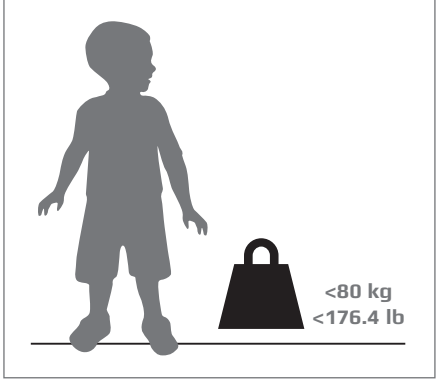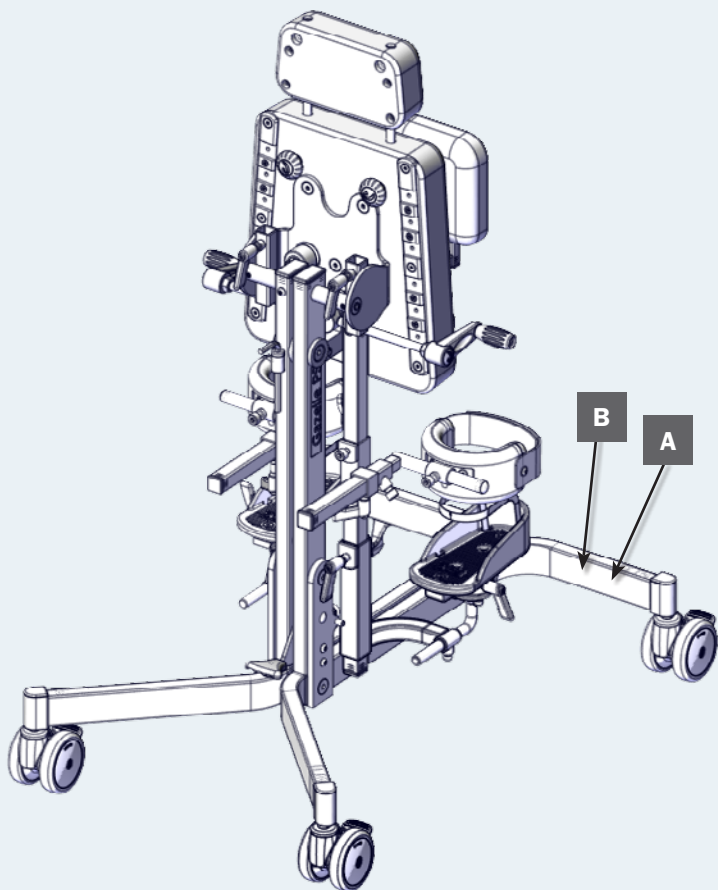A

R82[®]
by Etac

R82 A/S, Parallevej 3
DK-8751 Gedved
etac.com
Made in xxxxxxxx

MD

CE

UK
CA

YYYY-MM-DD

B

Posting group
Product name, size x

▲ Max. user mass: xxx kg
▲ Max. rated load: xxx kg
▲ Mass device: xxx kg




DMR ver.: xxxxx

SW xxxxxxxxxxxxxx

MP xxxxxxxxxxxxxxxxxxxx

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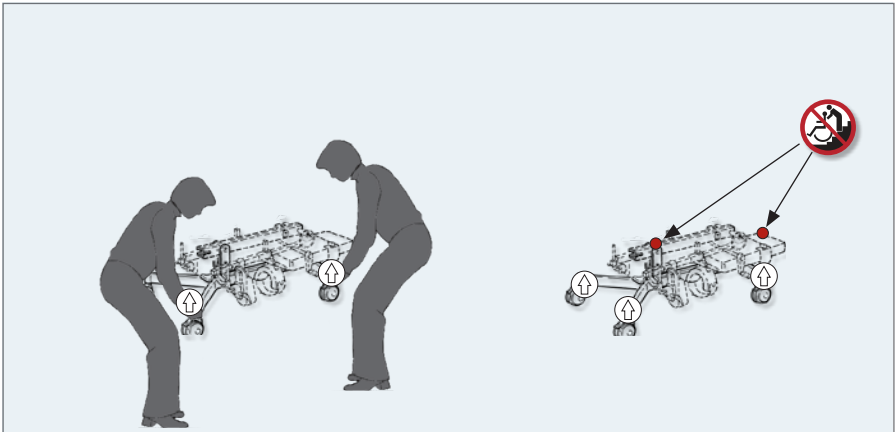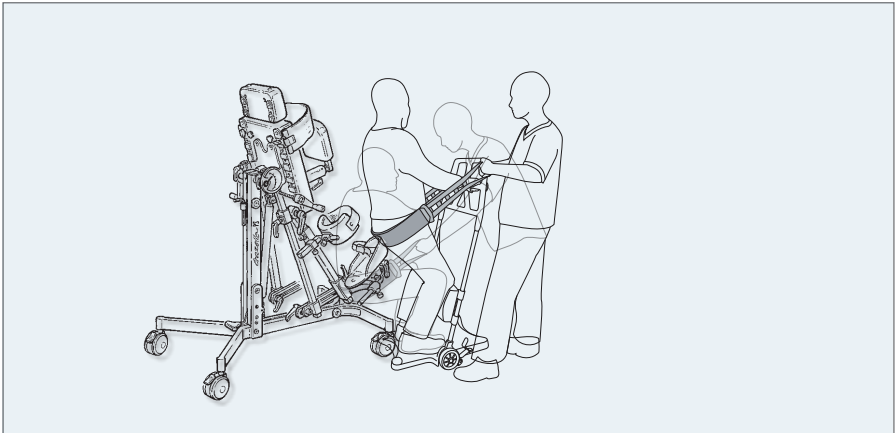
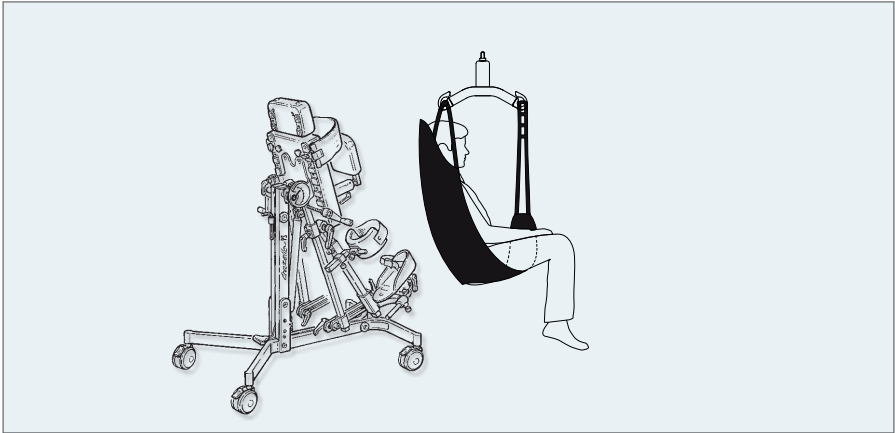
ENGLISH

Ready for use

The Gazelle PS is folded when delivered, to obtain as little space as possible, during transportation. Follow the description below, to unfold your Gazelle PS:

1. Make sure that the four wheels are locked.
2. Press against the floor frame with your foot, as the Gazelle PS is raised to vertical position. Make sure the Gazelle PS fits exactly into the fitting on the floor frame.
3. With the Gazelle PS in a vertical position, use the enclosed 6 mm Allen key to fasten the screw (A). The Allen key is placed on the Gazelle PS (B) and is always nearby for adjustments.

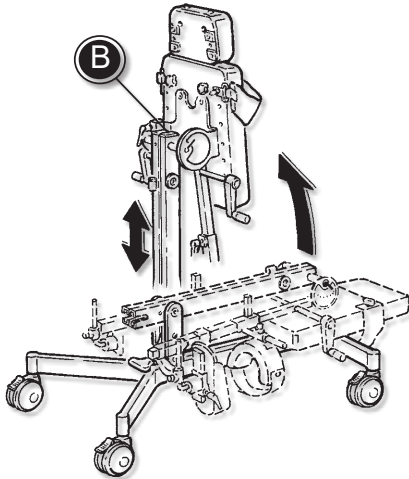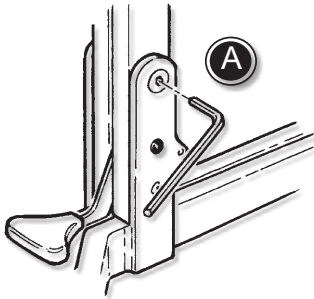**Head support with facial aperture**

Mounting of Head support with facial aperture:



- A. The two retaining bars go through the adjustable trunk support after removal of two plastic plugs.
- B. Positioning knobs can be screwed into place after making cut-outs in the adjustable trunk support.



It can be helpful to mount the head support with facial aperture before placing the child on the table in horizontal position. Remove the head support when the table is raised into vertical position.

Trunk support and extension

Height adjustable trunk support and extension can be mounted on the Gazelle PS (A).

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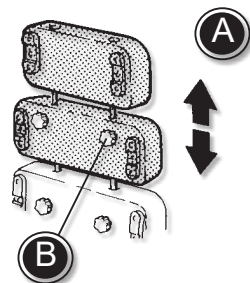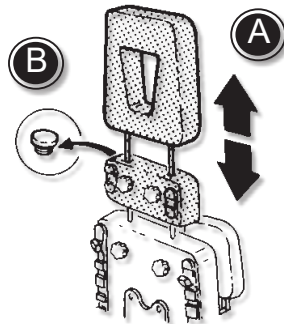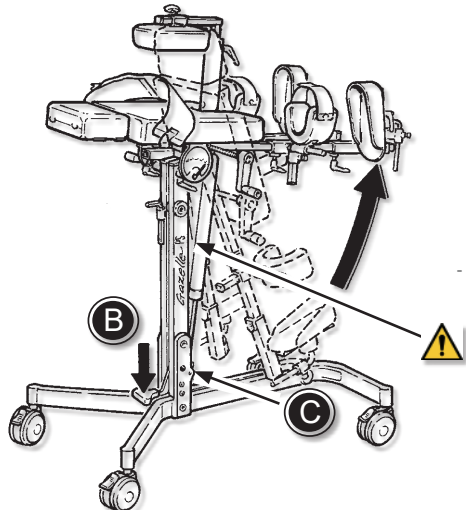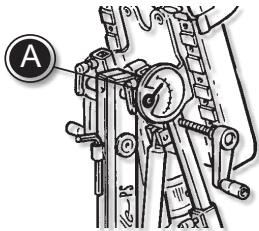Angling


Degree indicator (A)
Adjust the angle of the Gazelle PS from vertical to 15° under horizontal by pressing the foot pedal (B).
As security, the foot pedal can be locked at the wanted angle position. Turn the button at the lower end of the gasspring and the Gazelle PS locks (C).


- 05 Squeezing Hazard

Independent abduction

Adjusting the abduction 0° - 30° by turning the handle in both sides (A).
Retractable handle for tilt adjustment (B).

- 05 Squeezing hazard
- 05 Squeezing hazard

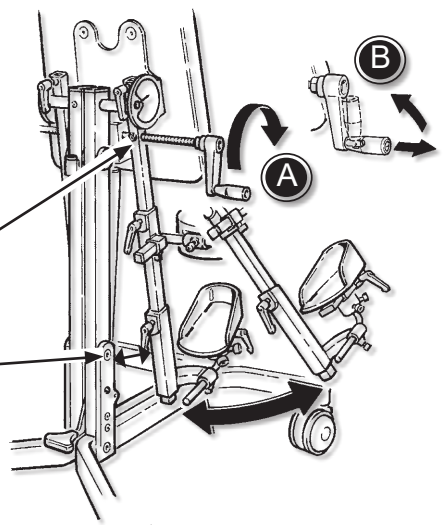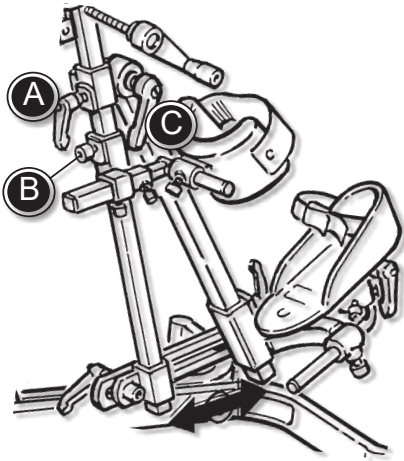Angling of lower leg

Mount the bar for the lower leg angling on the abduction bar, as follows:

1. Loosen the handle (A) to remove the support.
2. Loosen the Allen screw (B) to remove the knee supports.
3. Mount the knee support on the bar for the lower leg angling and tighten the handle (A).
4. Mount the foot support on the bar for the lower leg angling.


Adjust the height of the lower leg angling by turning the handle (B) and adjust the angling by turning the handle (C). To adjust the foot support see page 15.


Knee support

As standard, the Gazelle PS is delivered with knee supports mounted to be used by front-standing children. Turn around the knee supports, to be able to use it for backstanding children:

1. Loosen the handle (A) to remove the foot supports.
2. Loosen the handle (B) to remove the complete knee support unit.
3. Mount the complete knee support unit on the opposite side and fasten with the handle.
4. Turn the fitting for the knee support 180°, by loosening the screws (C) and (D).
5. Loosen the screw (E) and turn around the knee support, for backstanding children.

Adjust the knee support with the enclosed 6 mm Allen key. Loosen the screws (C) and (D) and the knee supports are depthadjustable. By loosening the screw (E) the knee supports adjust in width and angles. Use the handle (B) to adjust the height of the complete knee support unit.

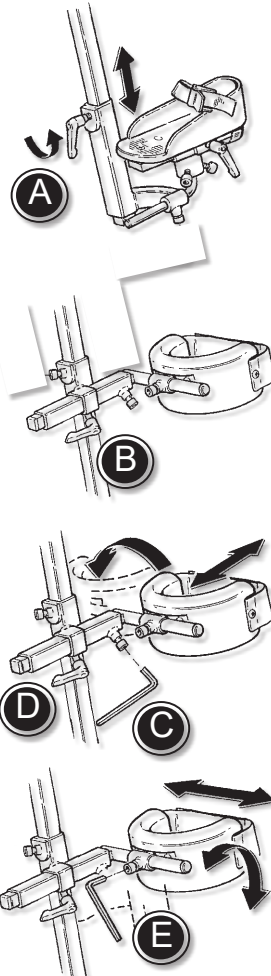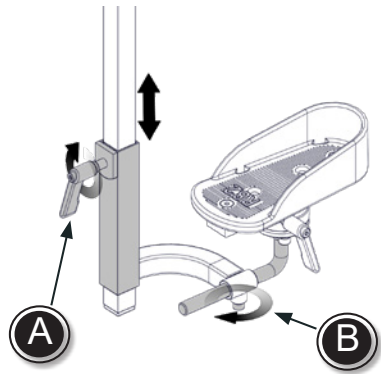Foot support





For use of back standing:

When the knee support is turned, the foot supports must be mounted in the same side as they were demounted.

1. Mount the foot supports and tighten the handle. (A).
2. Loosen the screw (B) and adjust the depth of the foot support.
3. Loosen the screws (C) and turn the foot support 180°.
4. Tighten the screws in the requested position.

The foot plate mounting is mounted on the joint. By means of this, the footplate can be adjusted in different angles by loosening the screws (C).

The plate can be adjusted in the depth by loosening the handle (D), or the unbracket screw (B).

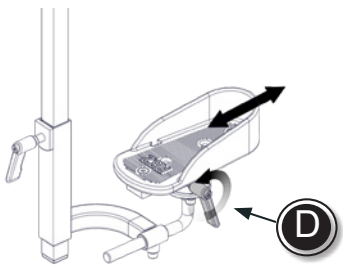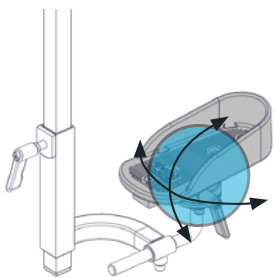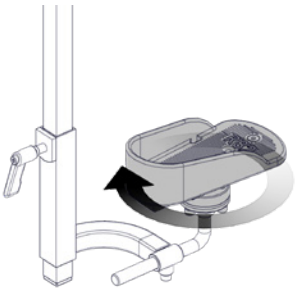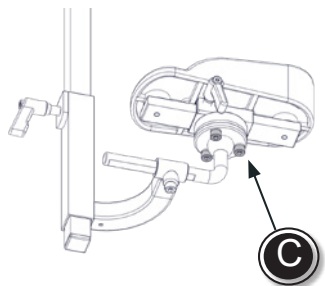
Strap

Mount the strap on the trunk support or on the basic plate considering where the support is needed.

The strap is mounted on the back (A) and fastened in front with velcro and a snap-lock fitting (B).

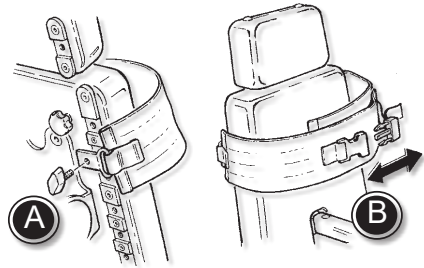Back support



The back support is adjustable in height (A) and in depth (B).

Before placing the child:

1. Remove the back support (C)
2. Fold down the fitting for back support by pulling the release chain (D).

Fasten with the handles.

The child may not sit on the mounting for the back support. There must always be a couple of centimetres between the child and the mounting.

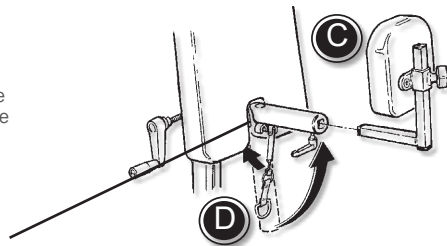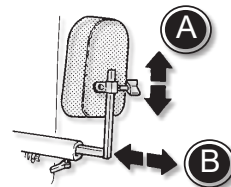05 Squeezing hazard!

Take care when refitting the back support.

Tray with / without bowl



Loosen the big handles (A) and (B) completely, and place the tray-mounting in the two tubes on the back. Fasten the big handles (A) and (B) again.

Loosen the big handle (A) and the tray-mounting can be turned through 180 degree's for either prone or supine use.

Place the tray in the tray-mounting tubes. When the required position is achieved, fasten the big handles (C)


Loosen the handles (B) to adjust the height of the tray. Loosen the handle (A) to adjust the angle for the tray.

The tray with bowl is delivered with an extra cover plate, so it can be used as a normal tray.

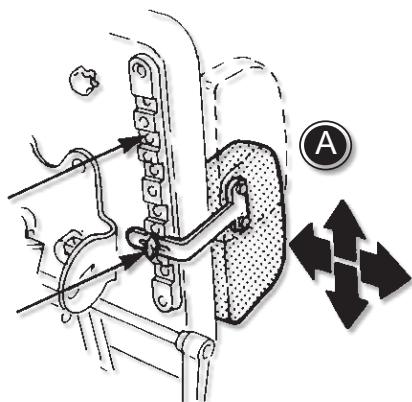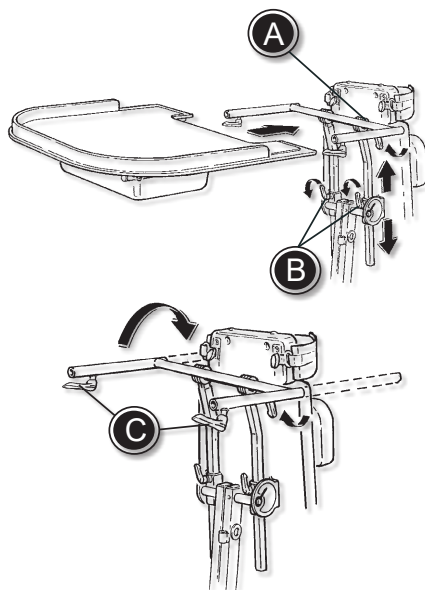Side support

Mount the side supports on the trunk support to give the rightsupport to the upper part of the body or on the basic plate to support the body and hips. The side supports are adjustable in height and in width (A).

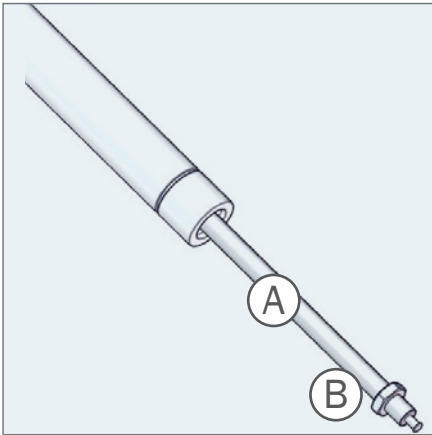
Adjusting the gas spring

Gazelle PS is equipped with a gas spring for angle adjustment. After a while it might need adjustment. The adjustment is necessary when:

- The gas spring is releasing itself
- The gas spring can not be released

If the gas spring releases itself, loosen the nut (B) with a 17 mm wrench. Turn the bar (A) anti-clock-wise to loosen. Tighten the nut (B) again.

If the gas spring will not release, loosen the nut (B) with a 17 mm wrench. Turn the bar (A) clockwise to tighten. Tighten the nut (B) again.

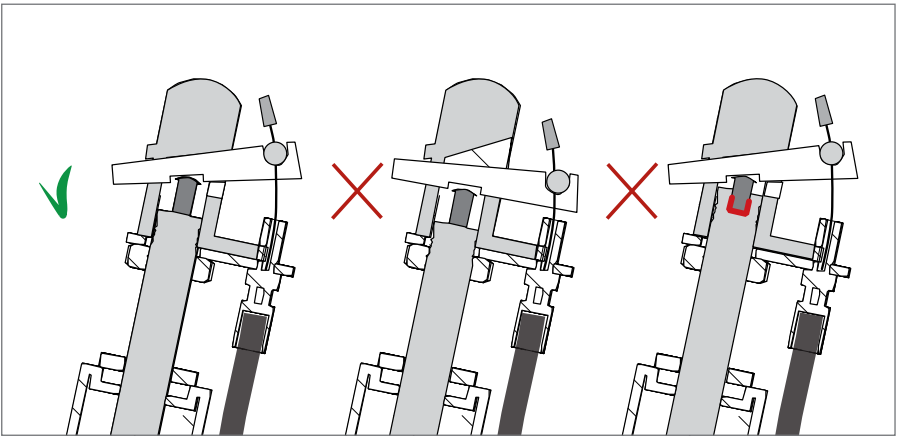General

We appreciate your choice of a new device from R82 - a worldwide supplier of technical aids and appliances for children and teenagers with special needs. To fully benefit from the options offered by this device please read this user manual before use and save it for future reference.

Intended purpose

The device is a standing frame for children and adults with disabilities.

Indications

The device can be used by children and adults with disabilities, who experience limitations in their ability to maintain a standing position. The target group for the device is based on functional ability and anthropometry. It is not based on a specific health condition or age. The user is the person standing in the frame. The caregiver is the person operating the frame. Never leave the user unattended in this product. Ensure permanent supervision by a caregiver. The use of the device will be discontinued more times during the day where the user will be lying down or sitting in another assistive device or a conventional product.

Contra indications

There are no known contra indications.

Precautions

Special attention should be given to the considerations about the configuration of the product and the use of accessories for users with limitation in their:

- Sensory functions and the sensation of pain, consider extra checks of the pressure caused by the support provided by the product and its accessories.
- Movement function (e.g. experience involuntary movements).
- Ability to maintain a standing position, use the product with accessories (e.g. chest or back support).
- Ability to maintain the position of their head, consider to use accessories or to use a product that provide support to the positioning of the head.

Application areas

The product is for indoor use

Enclosed tools:

- 6 mm Allen key

CE Declaration of Conformity

This device conforms to the requirements of the Medical Devices Regulation (2017/745). The CE mark must be removed when:

- Rebuilding the device
- Using the device outside the intended purpose and declaration of conformity
- Using other than original R82 spare parts and fittings.

Combinations

When combining the device with other devices, the combination must not change the intended purpose of the devices or modify the devices in such a way that compliance with the applicable requirements may be affected. The person / company responsible for the combination of the devices must ensure that the applicable requirements are fulfilled.

Applicable requirements examples:

- Material safety (biocompatibility, flammability), mechanical safety (stability, squeezing hazards, inter connection), electrical safety.

Disposal

When a product reaches the end of its lifetime, it should be separated by material type so that the parts can be recycled or disposed of properly. If required, please contact your local dealer for a precise description of each material. Please contact your local authorities to be informed of the options available for recycling in an environmentally-friendly way.

Accessories and spare parts

The products from R82 can be supplied with a variety of accessories including vests and belts which comply with the needs of the individual user. Spare parts can be ordered on request. Find the specific accessories on our website or contact your local dealer for further information.

R82 Warranty

www.etac.com/support

Presales information and instructions

Presales information and instructions are available at www.etac.com, under your device, in the documents section.

Safety

Signs, symbols and instructions placed on this product must never be covered or removed and must remain present and clearly legible throughout the entire lifetime of the product. Immediately replace or repair illegible or damaged signs, symbols and instructions. Contact your local dealer for instructions. In case of an adverse event occurred in relation to the device, incidents should be reported to your local dealer and the national competent authority in a timely manner. The local dealer will forward information to manufacturer.

Caregiver

- Read the "User Guide" carefully before use and keep it for future reference. Incorrect use of this product may cause serious injury to the user
- This product contains small parts which could present a choking hazard if removed from their intended placing
- Always use the correct lifting techniques and aids for this purpose
- Never leave the user unattended in this product. Ensure permanent supervision by a caregiver
- Repair/replacements must only be made using new original R82 spare parts and fittings and carried out in accordance with the guidelines and service intervals prescribed by the supplier
- If there is any doubt as to the continued safe use of your R82 product or if any parts should fail, stop using the product immediately and contact your local dealer as soon as possible
- Ensure that moving parts does not squeeze or trap body parts when adjusting them
- Ensure that body parts are not trapped in openings for mounting of supports/ accessories.
- Transportation of the product whilst the user is in it is not allowed. It is allowed to use the product with the user safely strapped in it for transfers over short distances. Do not cross any obstacles during transfer.

Environment

- Observe the surface temperature on the product, before placing an user in it. This goes especially for users with insensitive skin as they cannot feel the heat. If the surface temperature is above 41 °C then allow the product to cool before use
- Use this product on a level surface without any obstacles. Be aware of reduced maneuverability on slippery surfaces
- When using on inclines, please refer to measurements

User

- If the user is close to the max. load limit and/or has a lot of involuntary movements, e.g. rocking, a larger size with a higher max. load or a different R82

product should be considered

Product

- Carry out all positional adjustments on the product and accessories and ensure that all knobs, screws and buckles are securely fastened before use. Keep all tools out of the reach of children
- Activate the brakes before placing the user into the product
- Ensure the stability of the product is maintained before placing the user in it
- Ensure the castors are fully-functional and secure before use
- Inspect the product and all of its accessories and replace any worn out parts before use
- Do not expose the gas spring to pressure or high temperatures. DO NOT puncture
- Do not expose the piston rod to scratches etc. It is important to use a cloth or the like, if tools are used to adjust the gas spring
- The footplate is for support only! The foot support does not have a weight bearing function
- This product is tested for resistance to ignition at level non-ignition in accordance with EN 1021-1 and EN 1021-2
- The product might lose its flame resistant characteristics if using other than original R82 cushioning
- This product is not to be occupied in a motor vehicle. The user must be positioned in a vehicle seat and the product must be stored in a luggage compartment
- The user must not occupy the product while it is being transported into or out of a vehicle by a lift

Service & maintenance information

Service and maintenance of a medical device is wholly the responsibility of the owner of that device. Failure to service and maintain a device in accordance with the instructions may invalidate the device's warranty. Furthermore, failure to service and maintain a device may compromise the clinical condition or safety of users and/or their caregivers. Do not perform service and maintenance while the user is in the product. Contact the local dealer for assistance, if needed, in setting up, using or how to perform service and maintenance of the device.

Lifetime

The lifetime of this product in normal use is 8 years if all maintenance and servicing is carried out in accordance with the manufacturer's instructions and demonstrably recorded.

Reuse

This product and its accessories/components are suitable for reuse.

Refurbishment for reuse

Before reuse or changing the user, follow the service checklist, maintenance and washing instructions and carry out regular inspections.

Service interval

The frequency of inspection may be altered depending upon use and wear. It is recommended that the product is inspected yearly, every time the product is re-issued for use and after a long-term storage (more than 4 months). The inspection should be carried out by a person who understands the use of the product.

Service checklist:

Check and adjust the following as per their intended use:

- Fittings and buckles for signs of worn out parts.
- Castors are moving freely and all wheel locks/brakes are easily activated.
- Nuts and bolts on the product are tightened.
- Product for cracks or signs for worn-out parts.
- All labels on the product are intact.

Repair or replace any damaged or worn out parts.

Maintenance

Before every use it is recommended to: Wipe off any residue and dirt from the product, using a cloth with warm water and a mild detergent/soap without chlorine and let it dry before use.

Carry out a visual inspection, for any damaged or worn-out parts.

Washing

Main product

This product can be washed at 60° using a mild detergent for 10 minutes in a washing machine designed for medical devices. Use

the machine's drying function to dry the product. Unmount all accessories and wash them separately.

Handwash

This product and its accessories can also be washed by hand. Use warm water and a mild detergent/ soap without chlorine and let all parts dry before use.

Disinfection

The product can be disinfected with a 70% disinfectant IPA solution. It is recommended to wipe off any residue and dirt from the product, using a cloth with warm water and a mild detergent/soap without chlorine and let it dry before disinfection.

Materials

- Plastic
- Steel
- Fabric
- Wood


Surface treatment

The following surface treatment processes were used to protect against corrosion:

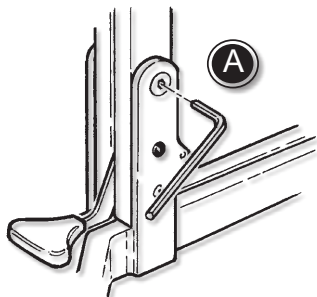- Lacquered surfaces = Polyester powder coating or ED-coating
- Non-lacquered aluminium parts = Anodized coating
- Non-lacquered steel surfaces = Galvanized

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Etac Ltd.
Gelbart AG